# **REF** TableTIM

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***BRUGSANVISNING***

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Product code: TableTIM



### 1. Characteristics

The bedside table is shaped like a C-shaped metal frame. It has an adjustable top on which objects can be placed and casters with a brake to ensure mobility and movement over the bed

#### Support components:

- a. Frame
- b. Base with wheels (2 sets)
- c. Tabletop frame
- d. Tabletop arbor
- e. Mounting butterfly
- f. Mounting butterfly with screw
- g. Bushing
- h. Screw with washer (2 sets)
- i. Tabletop
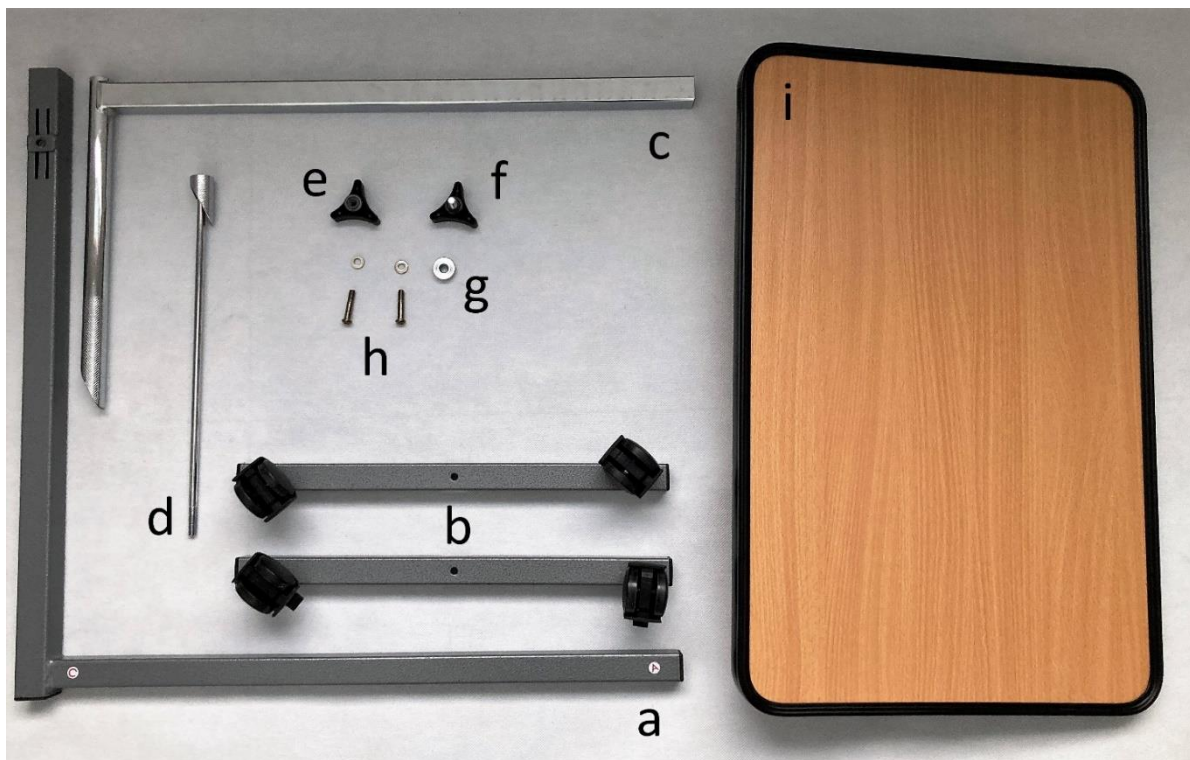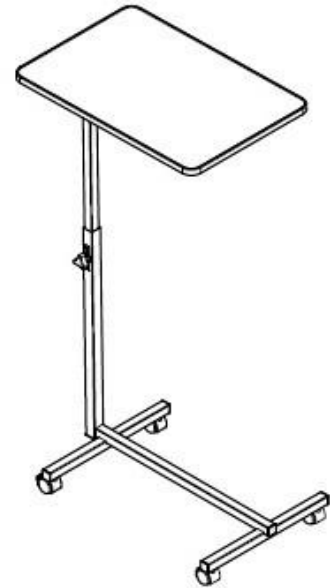
## **2. Application**

The bedside table is designed for use directly above a patient's bed to provide a surface for eating, reading or writing while in bed. The product is designed for use both during care in medical facilities and at home.

## **3. Assembly and adjustment**

- to the frame (a) screw both bases with wheels (b) - remember to properly match the letter markings - the letter "A" from the frame to the letter "A" of the base with wheels, similarly for the letters "B" / "C" / "D",
- after twisting the base, set it on 4 wheels and lock the wheel brakes for further assembly,
- into the square hole of the frame (a) insert the frame of the top (c),
- into the angled chamfered hole of the tabletop frame (c) insert the tabletop pin (d) in such a way that the two chamfered surfaces overlap,
- into the circular hole of the tabletop frame (c) insert the sleeve (g) on the protruding thread in such a way that the previously inserted tabletop mandrel (d) is not moved,
- on the protruding thread of the tabletop arbor (d) screw the mounting butterfly (e) - do not tighten the butterfly firmly - tighten it enough so that it does not fall out,
- put the tabletop (i) on the prepared tabletop frame. To facilitate the installation of the countertop, grasp it with one hand by the tube located from the underside of the countertop. After applying the countertop, tighten the mounting butterfly (e),
- fix the mounting butterfly with the screw (f) in the frame hole (a) and tighten it,

### **Adjusting the tilt of the tabletop**

- to adjust the angle of the tabletop, loosen the mounting butterfly (e) located on the side of the tabletop.
- after setting the correct angle, tighten the mounting butterfly (e)

### **Adjusting the height of the tabletop**

- to adjust the height of the tabletop, loosen the mounting butterfly with screw (f) located on the base frame (a).
- after setting the correct height, tighten the mounting butterfly with screw (f)

## Brakes

- lock the wheel brakes for the time of use of the table

## 4. Technical data

Dimensions (length x width)	60,5 cm x 40,2 cm
Tabletop dimensions	60 x 40 cm
angle adjustment of the tabletop	0° - 90°
Height adjustment	72 – 112 cm
Weight	6,7 kg
Max. capacity	10 kg
Wheels diameter	1,50"

## 5. Notes

**Timago International Group does not bear responsibility of improper use of the cane, failure to comply safety regulations and misuse.**

## 6. Maintenance and cleaning

The product should be cleaned with a damp cloth and then dried thoroughly. Do not use harsh polishing agents as this may damage the protective coatings. Contaminants should be regularly removed from the moving parts of the product (hinges, clamping screws, locks).

## 7. Storage and transport conditions

The product should not be exposed to direct sunlight, water or moisture.

## 8. Warranty information

All products distributed by our company are covered by a warranty, the terms of which are described in the warranty card available on our website. Please contact the retailer from whom you purchased the product. Please, keep in mind that for warranty purposes, the proof of purchase (receipt or invoice) should be retained.

## 9. Labels

	Reference number		Manufacturer
	LOT number		Manufacture date
	Serial number		Medical device
	Note		Please read the instruction
	The manufacturer checked for compliance with essential medical devices requirements.		

Our complete line of products is available from our distributors.

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